



iM70
Patient Monitor

iM70 Patient Monitor



About Edan

Edan is a healthcare company dedicated to improving the human condition around the world by delivering value-driven, innovative and high-quality medical products and services. For over 20 years, Edan has been pioneering a comprehensive line of medical solutions that address a broad range of healthcare practices including:

- Diagnostic ECG
- Patient Monitoring
- OB/GYN
- Ultrasound Imaging
- Point-of-Care Testing
- In-Vitro Diagnostics
- Veterinary

Healthcare professionals around the world depend on Edan's breakthrough medical technologies and outstanding customer support.

EDANUSA

4901 Morena Blvd., Suite 505 | San Deigo, CA 92117 | Tel: 888.850.4597

EDANUSA
www.edanusa.com

iM70

Patient Monitor



12.1"

Maximum 11 Waveforms




No-fan Design Night Mode



Pitch Tone Thermal Recorder



Shortcut Menu

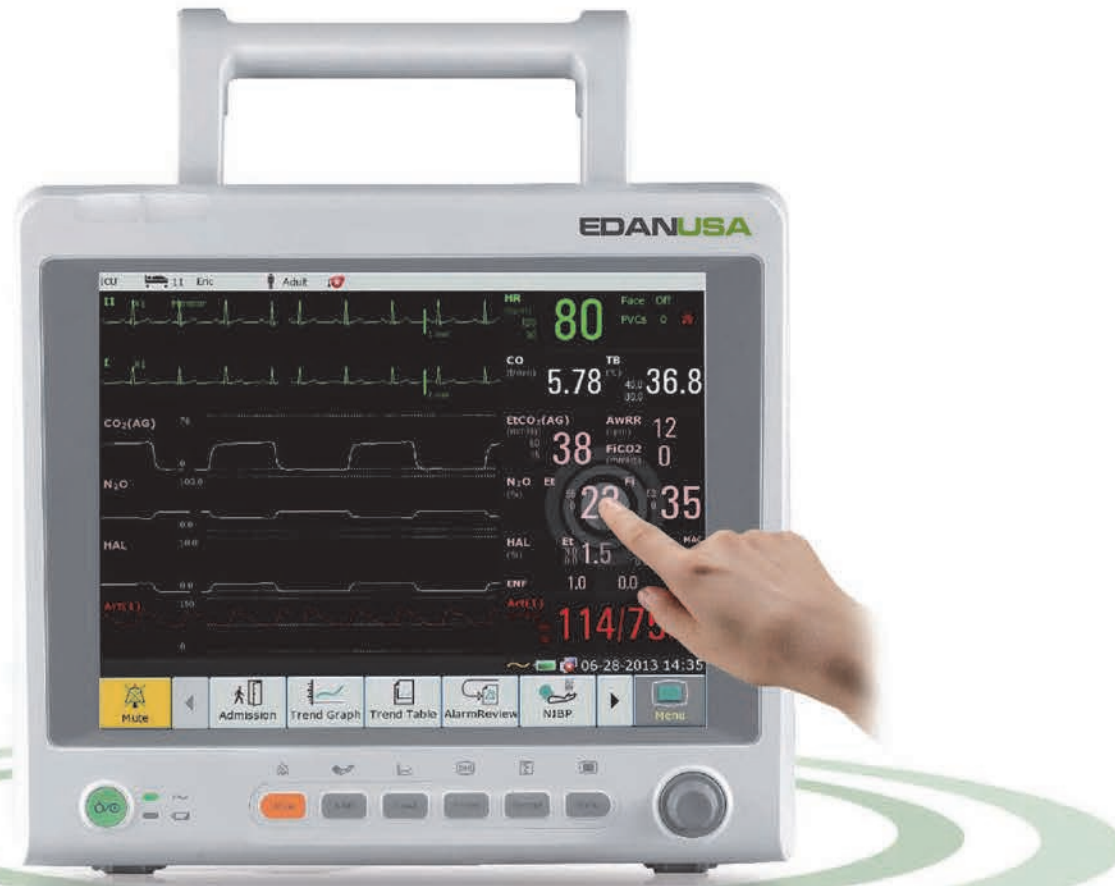

Standard Parameters

- 3/5-lead ECG, HR, RESP, EDAN SpO₂, EDAN NIBP, PR, 2-TEMP

Optional Parameters

- 2-IBP, C.O., EDAN G2 CO₂ (Sidestream), Respironics CO₂ (Mainstream and Sidestream), MASIMO Anesthesia Gas/O₂, Nellcor OxiMax™ SpO₂, Omron® NIBP



Reliable Algorithms

- iSEAP™ ECG algorithm optimized for arrhythmia detection, pacemaker detection, and HR measurement
- iMAT™ SpO₂ algorithm with outstanding motion resistance and low perfusion resistance performance
- ICUFS™ NIBP algorithm optimized for cardiac patients, hypertensive patients, and neonatal patients

Multiple Display Modes



End tidal Carbon Dioxide for Intubated/Non-Intubated Patients



EDAN G2 CO₂ (Sidestream)

- Superior water trap design for accurate monitoring
- iCARB™ algorithm with intelligent CO₂ pseudo wave identification technology
- Multiple sampling accessories as options for adult, child and neonate patients

Respironics CO₂ (Mainstream/Sidestream)

- Plug & play module design
- Dehumidification tube instead of water trap
- Low sampling rate of 50ml/min suitable for all types of patients

MASIMO Anesthetic Gas/O₂ (Mainstream/Sidestream)



- Unique Plug & Play module for mainstream anesthetic gas monitoring
- The no moisture sampling line design suitable for all clinical applications
- A stable low flow of 50ml/min for all types of patients



Nurse Call



VGA Output



Defib Proof & Sync



ESU Proof



Bed Rail Hook



Wall Mount



Rolling Stand

Built-in Temporary Memory:

1200 NIBP Measurement Review	120 h Trend Review	60 Alarm Review	120 s Frozen Waveform
---------------------------------	-----------------------	--------------------	--------------------------